



L&P MEDICAL TEL.0740524960	VERENDEANU IONEL EMAIL:verendeanuionel@gmail.com	
RO42026251		

KIT CORONARICO/CORONARY SET

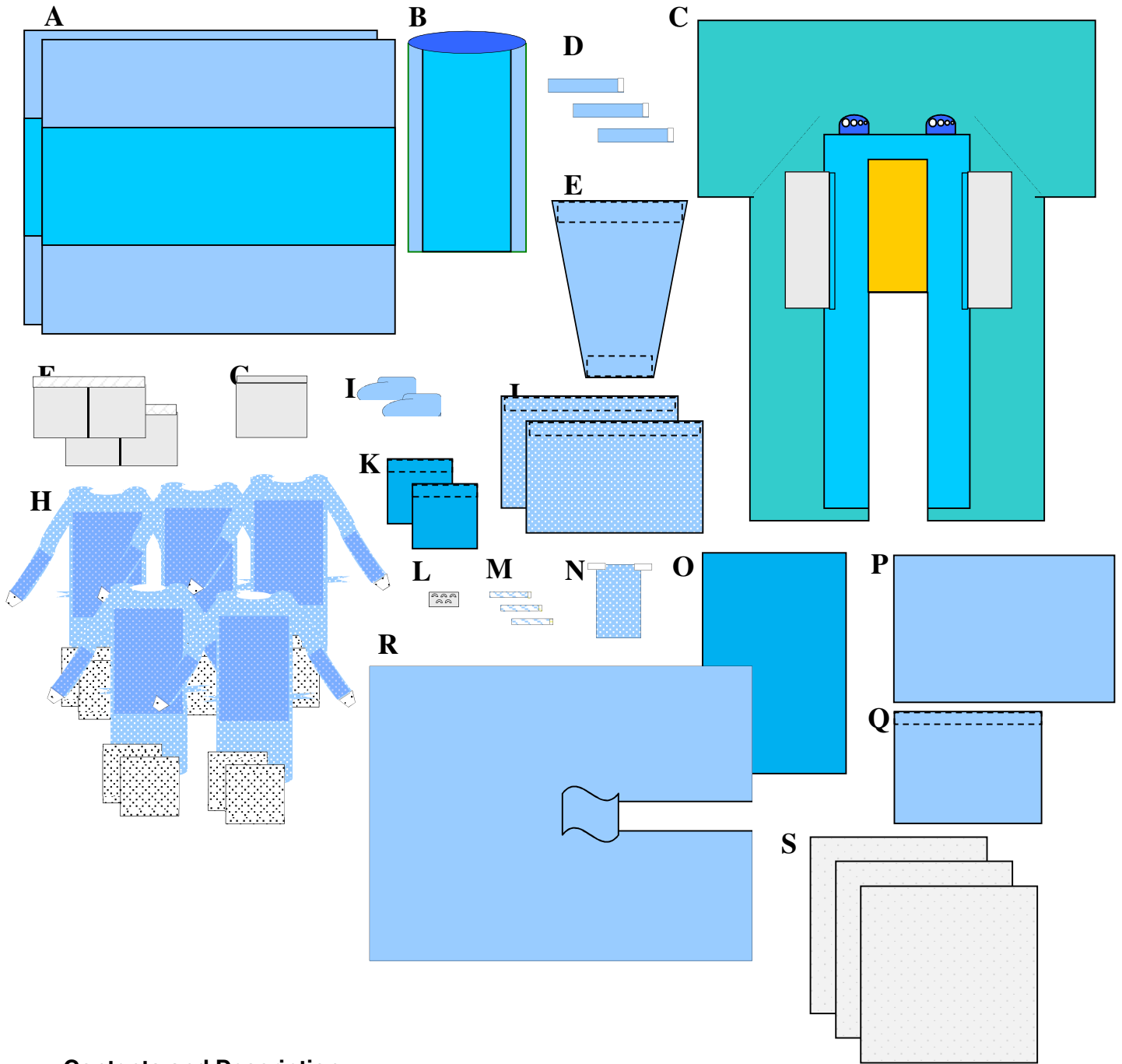

Contents and Description:

- A) 2 200x240 cm instrument cover table drape, one of this wrapping (32.600)
- B) 1 80x145 cm mayo stand cover (32.002)
- C) 1 270/200x335 cm "T" patient drape with 41x82 cm "U" split aperture covered by adhesive incise film with iodine solution, super absorbent reinforcement, lateral instrument pouch, multi holes guide for tubes and cables passage (32.141-I)
- D) 3 10x50 cm adhesive strips with double finger opening (33.023)
- E) 1 75x150 cm mayo/operating table cover drape with adhesive (31.120/SRV)

Data 20/10/20	COMPILED BY: 	Pagina 1 di 3
---------------	---	---------------

	TECHNICAL FILE REF. MTT0635-I/02	M-750-02 REV.C DEL 18/09/19
		

- F) 2 35x70 cm double compartments instrument pouch with adhesive (33.054/BET)
- G) 1 35x40 cm waste bag with fixing adhesive (33.019)
- H) 5 HSB wrapping surgery gowns size XXL with waterproof reinforcement In the thorax and in the sleeves, heat sealed sleeves, Velcro neck closure for security, 4 waist belts for inside and outside closure, turn-up folding, planned opening with 2 30x40 cm absorbent towels (33.046)
- I) 2 shoe covers with elastic (A951N/M)
- J) 2 75x120 cm lateral cover drapes with long adhesive side (31.229)
- K) 2 45x45 cm super absorbent drapes with fixing adhesive (31.097/SRV)
- L) 1 tubing holder with 5 cable ties and fixing adhesive (33.135/1)
- M) 3 2.5x12 cm velcro tapes with fixing adhesive (15.003/TR)
- N) 1 30x50 cm shaped pubis cover drape with fixing adhesive (33.150/SRV)
- O) 1 97.5x150 cm reinforcement instrument cover drape (14.019/ES)
- P) 1 100x150 cm drape (31.107)
- Q) 1 75x100 drape with long adhesive side (31.120)
- R) 1 200x260 cm with 20x100 cm "U" split adhesive aperture with panties (32.026/K-4)
- S) 3 120x120 cm anti-contamination wrapping envelopes (13.058)

Use Destination:

- Sterile patient protection and surgical field isolation during cardio surgery procedures

Features:

- Three wrapping envelopes and anti-contamination protective table drape
- Removable labels for traceability
- Latex free
- Suitable bending to guarantee the total respect of the corrected aseptic technical in the creation of the operating field

Materials used:

- Table drape and mayo made of polyethylene 45 g/m² and nonwoven reinforcement 35 g/m² (A,B)
- Patient drape made of hydro-repellent azure nonwoven spunlaced 70 g/m² with super absorbent reinforcement 110 g/m² incise film in antistatic medical PU 25 µm (C)
- Strips made of laminated adhesive nonwoven 60 g/m² (D)
- Drape made of azure two layers nonwoven 53 g/m² with medical adhesive polypropylene 50 g/m² acrylic base, biocompatible (E,Q,P)
- Pouch made of antistatic polyethylene 60 µm (G,F)
- Gown made of hydro-repellent nonwoven azure SMMS 50 g/m² with waterproof reinforcement and towels made of absorbent cellulose 50 g/m² (H)
- Shoe cover made of azure CPE 12 g/m² with latex free elastic (I)
- Drape made of hydro-repellent azure nonwoven SMS 50 g/m² with medical adhesive polypropylene 50 g/m² acrylic base, biocompatible (J,N)
- Drape made of super absorbent nonwoven 113 g/m² with medical adhesive polypropylene 50 g/m² acrylic base, biocompatible (K)
- Tubing holder made of high density sponge, hook in reinforced PA (L)
- Strips made of white Velcro with medical adhesive (M)
- Drape made of super absorbent nonwoven 113 g/m² (O)
- Drapes made of azure three layers nonwoven 73 g/m² with medical adhesive polypropylene 50 g/m² acrylic base, biocompatible (R)
- Envelopes made of white crepe paper 52 g/m² (S)
- There are not incompatibility of the above written materials with fabrics and biological liquids
- Specific Data Sheets of medical devices are available on request

Data 20/10/20	COMPILED BY: 	Pagina 2 di 3
---------------	--	---------------

	TECHNICAL FILE REF. MTT0635-I/02	M-750-02 REV.C DEL 18/09/19
		

Storage and dispose:

- The device should be stored in cool place away from sources of heat
- The product is completely disposable like refusal hospital

Sterilization:

- Ethylene Oxide. Shelf life: 3 years

Conformity:

- CE marking, 93/42/CEE directive. Notified Body: 0373
- Device manufactured according to internal procedures of the ISO 13485 and ISO 9001 certified quality management system
- Woven Non woven validated according to EN13795 1-2-3 harmonized standards
- Class sterile III
- CND C99-GMDN 46370
- This device is custom made. The registration number will be provided at the market introduction time
- Sterile packaging and shelf life validated according to E11607 1-2 harmonized standards

Packaging:

- Single packaging in sterile TYVEK/PE bag 61x76 cm
- Double package inside shipping carton for safe near-O.R. storage
- Shipping carton 59x39x33 cm contains 1 piece

Data 20/10/20	COMPILED BY: 	Pagina 3 di 3
---------------	--	----------------------